

LIMITED PRODUCT WARRANTY

Newline warrants its products to be free from defects in material and workmanship for a period of one year from the date of shipment. Any product found to be defective within this time period may be returned to Newline, freight prepaid by buyer, with prior return authorization for repair or replacement at no charge. Any damage on the products which resulted from mishandling, physical abuse, accident, unauthorized service or product abuse will have the warranty voided immediately.

Newline's liability under this warranty is limited to the repair or replacement of the defective product and in no event shall Newline be liable for consequential or indirect damages to equipment or personnel. Nor shall Newline be liable for damages to equipment or personal injury caused by misuse, overload, accidental damage, alteration, improper installation or unauthorized opening of the equipment. Under no circumstances will Newline be responsible for any indirect or consequential damages due to errors in weighing or failure of a Newline product to perform properly.

To request warranty service, please email your claim with invoice # to our Service Department, service@newlinescale.com, to request a RMA#. You can also fill out the RMA form from our web site www.newlinescale.com



Super Mini Bathroom Scale

INSTRUCTION

For the purpose of proper using & maintaining the scale, please read this instruction completely before operating.



Super Mini Bathroom Scale

Your NewlineNY Super Mini Travel Bathroom Scale has been designed to accurately indicate your weight gain or loss over a period of time, and should give many years of service under normal use.

SPECIFICATIONS

- Equipped with a high precision strain gauge sensor system
- Capacity: 150 kg / 330 lb
- Division: 0.1 kg / 0.2 lb
- Automatic zero resetting / switch off
- Unit switch: kg / lb / st
- Low battery / over-load indicator
- Newly improved Auto Step-On Technology (no need to foot-tap activation before stepping -on)

POWER SUPPLY

- Battery: 1 x CR2032

OPERATION

STEP 1: For lithium battery powered scale, please pull the insulation tape out before first use.




STEP 2: When the LCD indicate "Lo", open the battery door and replace the battery with a new one.



STEP 3: This scale must be operated on flat, level & solid floor for accurate results. Do not use it on carpet or rugs.

STEP 4: Step on the scale evenly without moving, the scale will automatically display accurate weight in KG, LB or ST instantly. (No need to foot-tap to turn on the scale).

STEP 5: Put your feet evenly on the scale without moving; when the indication light is red, it means unstable. If you stand firm, the indication light will be green. Then, you can step off from the scale, your correct weight should be displayed.

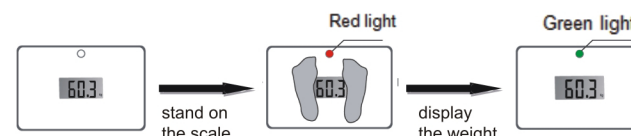


Put the scale on flat & level floor Start to measure the weight Finish the weighing

AUTOMATIC SWITCH OFF

The scale will be switched off automatically after 6 seconds of idle time.

UNIT SWITCH

The scale by default is set to Kg. Please press the switch button on the back of the scale to switch to Lb or St.

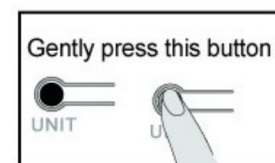

OVER-LOAD INDICATOR

When the scale is over 150 kgs, the "Err" will be displayed. Frequently over weight will also damage the sensor.



ADVICE FOR USE AND CARE

- Place your scale on a flat, level & hard surface will ensure accuracy of your scale.
- Always weigh & place your scale on a same spot.
- Do not let water get into the scale, as it will damage the electronics
- Do not use this scale for legal trade
- Treat the scale with care, as it is a precision instrument. Do not drop or jump on the scale to avoid damaging the sensor