



EXCLUSIVE
SensorDisc
TECHNOLOGY

LONG LIFE
BATTERY

RoHS

NEWLINE[®]

Electronic Bathroom Scale

INSTRUCTION

Item No. **SBB0807**

For the purpose of proper
using & maintaining the scale,
please read this instruction
completely before operating.

Electronic Bathroom Scale


Your new electronic personal scale has been designed to accurately indicate your weight gain or loss over a period of time and should give many years of service under normal use .

SPECIFICATIONS

- Equipped with a high precision strain gauge sensor system
- Capacity: 180 kgs/400lb
- Division: 0.1kg/0.2lb
- Automatic zero resetting / switch off
- Unit switch: kg / lb / st
- Low battery / over-load indicator
- Newly improved Auto Weigh & Display function (no need to foot-tap activation before stepping -on)

POWER SUPPLY


- Battery: 1xR2032
- For lithium battery powered scale, please pull the insulation tape out before first using.
- When the LCD indicate "Lo", open the battery door and replace the battery with a new one.

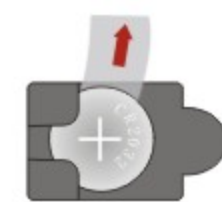

OPERATION

STEP1: This scale must be operated on flat & solid floor for accurate results. Do not use it on carpet or rugs.

STEP2: Step on the scale evenly without moving, the scale will automatically display accurate weight in KG, LB or ST instantly. (No need to foot-tap to turn on the scale).

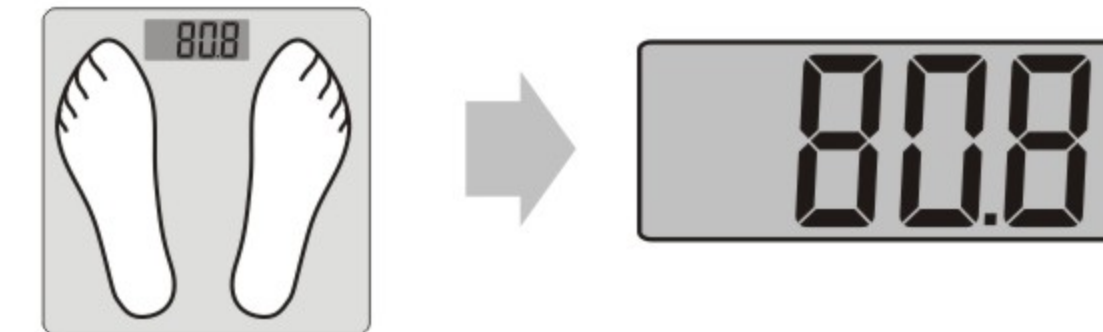

AUTOMATIC SWITCH OFF

The scale will be switched off automatically after 6 seconds of idle time.

OVER-LOAD INDICATOR

When the scale is over loaded(over 180kg) "Err" will appear on the screen:



ADVICE FOR USE AND CARE

- Place your scale on a solid & flat surface will ensure accuracy of your scale.
- Always weigh & place your scale on a same spot
- Do not let water get into the scale, as it will damage the electronics
- Do not use this scale for legal trade
- Treat the scale with care, as it is a precision instrument. Do not drop or jump on the scale to damage the sensor